Sapphire

Anterior Cervical Plate System



Confidence With Every Turn.



Sapphire® | Anterior Cervical Plate System

- Low profile, sleek, and narrow plating system for optimal anatomical fit
- Provides constrained, semi-constrained, or hybrid screwconfigurations
- Large plate windows allow exceptional graft visualization
- Integrated locking system for confident results
- Variety of preparations and screw placement options
- Reliable and intuitive instruments for simplified implantation



Thin Plate Profile

Sapphire is one of the lowest profile plating systems on the market. The Sapphire plate is 2mm thin, minimizing protrusion.

Contoured Design

The Sapphire plate is prelordosed and narrow to fit with most anatomical curvatures. It offers surgeons the flexibility to further contour the plate at the "bend zones" for unique patient anatomy.

Large Windows

Sapphire's windows allow for exceptional visualization of the interbody space.





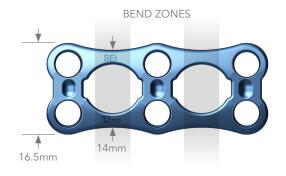


Plate System Angulation

- 7° Built into the plate
- 8° Additional screw angulation (+/-)
 - 15° Maximum angulation from center of plate curvature

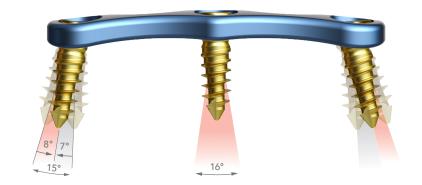

Sapphire® | Anterior Cervical Plate System

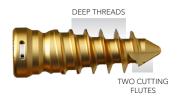

UNLOCKED NITINOL CLIP LOCKING SCREW BONE SCREW LOCKING MECHANISM

Proven Screw Design

Sapphire's screw design allows for excellent screw retention, robust thread engagement, and easy insertion.

Integrated Locking System

Sapphire's locking mechanism is built into the head of the bone screw. It is intended to help resist screw backout, maintain the low profile of the plate, and simplify the surgical technique.

Audible and Tactile Locking Confirmation

Once the locking mechanism has been deployed an audible click can be heard and felt, delivering confidence.


Intuitive Instrumentation

The Sapphire system can simplify procedures with its functionally elegant instrumentation.

- The Screwdriver design provides a secure interface with screws.
- The Drill Guides lock into plates for positive fit.



Screwdriver Unique four-lobe design for secure retention.

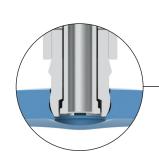

Drill Guides
Lock into plate
automatically.



Sapphire® | Anterior Cervical Plate System





Sapphire® Plate

Level	Length
1	11–25
2	26-44
3	42-68
4	60-84
5	85-105

Sapphire® Screws

Sappinic Sciews		
Diameter 4.0 4.5	Length 10 Magenta 12 Gold 14 Blue 16 Green 18 Bronze	Type Variable: Self-Tapping Self-Drilling* Fixed: Self-Tapping Self-Drilling*

Dimensions expressed in millimeters.
* Self-Drilling available in 4.0 only and 10–16mm.

